- Ⓛ DE Bedienungsanleitung
- Ⓛ GB Operating Instructions
- Ⓛ FR Mode d'emploi
- Ⓛ NL Handleiding
- Ⓛ IT Istruzioni per l'uso
- Ⓛ ES Instrucciones de uso
- Ⓛ PT Manual de utilização
- Ⓛ RU Руководство по эксплуатации

DE	Bedienungsanleitung	4
GB	Operating Instructions	7
FR	Mode d'emploi	9
NL	Handleiding	12
IT	Istruzioni per l'uso	15
ES	Instrucciones de uso	18
PT	Manual de utilização	21
RU	Руководство по эксплуатации	24

1 2 3 4 5 6

7 8 9 10 11 12

13

14

15

16

18

19


20

21

22

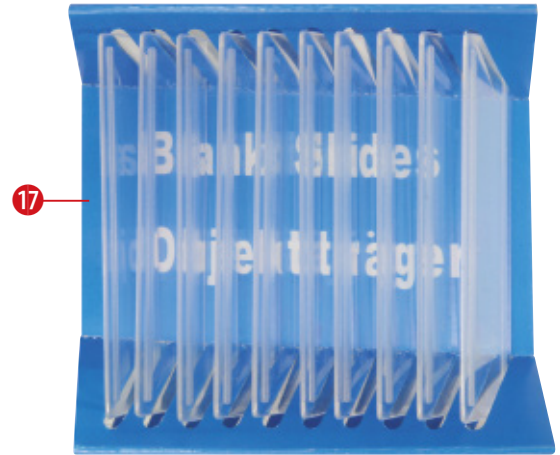17

Blank Slides
Objektträger

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RISK to your child!

Aids with sharp edges and tips are sometimes used with this device. Please store the device and all of its accessories and aids out of the reach of children. There is a risk of **INJURY**.

Children should only use this device under supervision. Keep packaging materials (plastic bags, rubber bands, etc.) away from children. There is a risk of **SUFFOCATION**.

The chemicals and liquids provided should be kept out of reach of children. Do not drink the chemicals! Hands should be washed thoroughly under running water after use. In case of accidental contact with the eyes or mouth rinse with water. Seek medical treatment for ailments arising from contact with the chemical substances and take the chemicals with you to the doctor.

Fire/Burning RISK!

Never subject the device - especially the lenses - to direct sunlight. Light ray concentration can cause fires and/or burns.

TIPS on cleaning

Clean the lens (objective and eyepiece) only with the cloth supplied or some other soft lint-free cloth (e.g.

micro-fibre). Do not use excessive pressure - this may scratch the lens.

Dampen the cleaning cloth with a spectacle cleaning fluid and use it on very dirty lenses.

DISPOSAL

Dispose of the packaging material/s as legally required. Consult the local authority on the matter if necessary.

Your experiment set is made up of the following components:

- 1 Shrimp eggs
- 2 Sea salt
- 3 Yeast
- 4 Gum media
- 5 Prepared slide "Fabric" (material)
- 6 Prepared slide "Shrimp"
- 7 Prepared slide "Honeybee"
- 8-12 Empty vials
- 13 Prepared slide "Butterfly Wings"
- 14 Cover slips
- 15 Adhesive labels
- 16 Shrimp hatchery
- 17 Slide
- 18 Pipette
- 19 Tweezers
- 20 Dissecting needle
- 21 Scalpel
- 22 Magnifying glass

What you need to know

Caution!

When using your experiment set, wear an apron (or old clothing) and rubber gloves!

Always carry out your experiments under the guidance and with the help of an adult!

Your experiment set includes numerous accessories that help you carry out different experiments. You can use all the parts independently of the microscope that you use for your research. Use the following web link to find interesting experiments you can try out.

<http://www.bresser.de/c/de/downloads>

You can do them with the individual parts of your experiment set.

Shrimp eggs

These are special eggs that are capable of survival even when they are dried out. As a result, they are particularly ideal for raising brine shrimp. You can find out how in the appendix to these instructions.

Sea salt

You can observe the sea salt using your microscope. Salt is a crystal that looks very interesting through a microscope. Apart from that, you will also need the sea salt to make a

salt solution, in which you will hatch your brine shrimp.

Yeast

The yeast in your experiment set is included as feed for the brine shrimp. If you feed them regularly, you will be able to watch how they grow.

Gum media

This material is a glue that you can use to make permanent prepared specimens.

Prepared slides „Honeybee,” “Fabric” (material), “Shrimp” and “Butterfly Wings”

Different prepared slides that you can investigate using your microscope.

Empty vials

You can keep a few samples in the empty vials (e.g. leaves, and, etc.)

Cover slips

You use the cover slips to cover specimens that you have placed on a slide in order to observe with a microscope.

Adhesive labels

You can use these labels to identify your permanent prepared specimens. You can also use them to label the empty vials.

Shrimp hatchery

You can hatch shrimp in the different chambers and observe the different stages of development.

Slides

You place the specimen on a slide. After you have covered it with some water or gum media, you can place a cover slip on top. Then you place the slide under the clips on the microscope.

Scalpel

The scalpel is a sharp knife that you can use to slice specimens/samples.

Tweezers

These are a type of pincers that you can use to pick up smaller objects.

Dissecting needle

A dissecting needle can be used for many things. You can use it to scratch the surface of a specimen, to fix objects in place or to stir liquids.

Pipette

You can use this instrument to soak up liquids in small amounts and squirt them out again.

To soak up liquids, push the upper part (head) of the pipette together and put the opening on the lower part in the liquid. Let the head go and the pipette will fill with the liquid.

When you want to release some of the liquid, simply push the head together again. The more carefully you press, the less liquid will come out of the opening.

Magnifying glass

With the magnifying glass, you can observe objects with 2 times the magnification.

After use

To make your experiment set - and the enjoyment of it - last longer, you should make note of the following for after each use.

1. Clean the cover slips, slides and the microscope instruments carefully with water and some soap.

Caution!

Cover slips, slides and the instruments have sharp edges and are in some cases pointed. Make sure to always let adults help you when cleaning, so that you don't hurt yourself!

Warning! The cover slips are particularly thin and could break. It is best to place them on a piece of paper on a hard surface for cleaning.

2. Close all the vials tightly after use. Some containers have liquids that evaporate very easily, and could dry up.

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Bresser GmbH